



Q2ROLLER®

Lateral Turning Device

User Manual




















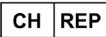


Visit www.HoverMatt.com for other languages

TABLE OF CONTENTS

Symbol Reference	2
Intended Use and Precautions.....	3
Part Identification - Q2Roller	3
Part Identification - Air Supply	3
Air Supply Keypad Functions.....	4
Instructions for Use	4
Product Specifications/Required Accessories.....	5
Cleaning & Maintenance	6
Returns and Repairs.....	7

Symbol Reference

	CE MARKING OF CONFORMITY		LATEX FREE
	UK MARKING OF CONFORMITY		LOT NUMBER
	AUTHORIZED REPRESENTATIVE		MANUFACTURER
	UK responsible person		DATE OF MANUFACTURE
	Switzerland AUTHORIZED REPRESENTATIVE		MEDICAL DEVICE
	CAUTION / WARNING		SINGLE PATIENT - MULTIPLE USE
	IMPORTER		Do not launder
	DISPOSAL		UNIQUE DEVICE IDENTIFIER
	OPERATING INSTRUCTIONS		PATIENT WEIGHT LIMIT
	USE TWO CAREGIVERS		
	LOCK ALL WHEELS		

Intended Use and Precautions

INDICATIONS FOR USE

- Assist patients who are unable to turn or laterally position themselves.
- Patients that need Q2 turning for off-loading pressure.
- Use as a Post-Prone Positioning Device in lieu of pillows, foam, or gel positioners.
- Assist patients to sit on the edge of the bed.

CONTRAINDICATIONS

Patients who are experiencing thoracic, cervical or lumbar fractures that are deemed unstable.

INTENDED CARE SETTINGS

Hospitals, long-term or extended care facilities.

PRECAUTIONS

- Ensure the Q2Roller meets the patient's needs and use clinical judgement in using the Q2Roller.
- It is recommended to use a minimum of two caregivers when operating the Q2Roller. However, one caregiver may operate the Q2Roller if side rails are raised.
- Ensure valves are not in contact with the patient while on the Q2Roller.
- Product is single-patient use and should only be used by trained personnel.

- Only use attachments and/or accessories that are authorized by HoverTech.
- Use the Q2Roller according to instructions.
- Make sure patient is centered on the Q2Roller before inflating using the Reference Center Line Indicator on the Q2Roller.
- Caregivers must ensure patient is attended and stabilized during Q2Roller inflation, deflation, or when inflated at angles greater than 30 degrees.
- Do not launder.
- Reference product-specific user manuals for additional operating instructions.



Only for use with the HT-Air air supply adjustable button.


PRECAUTIONS – HOVERTECH AIR SUPPLY

- Route the power cord in a manner to ensure freedom from hazard.
- Avoid blocking the air intakes of the air supply.
- Use this product only for its intended purpose as described in this manual.




Avoid electric shock. Do not open air supply.

Part Identification - Q2Roller

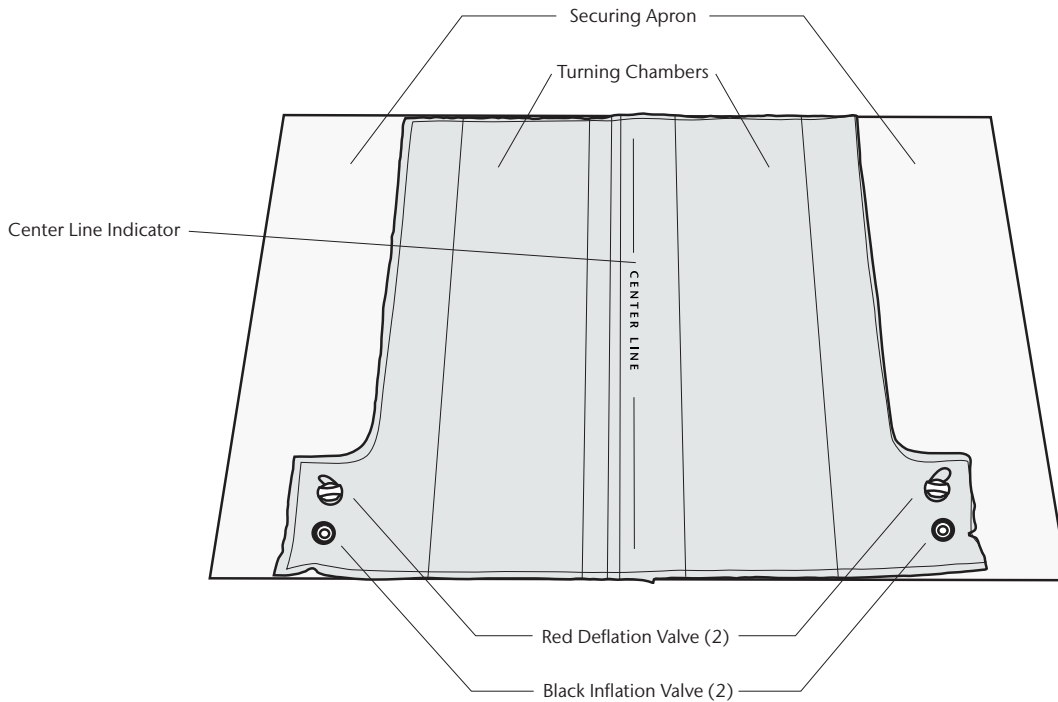


Part Identification – HT-Air® Air Supply

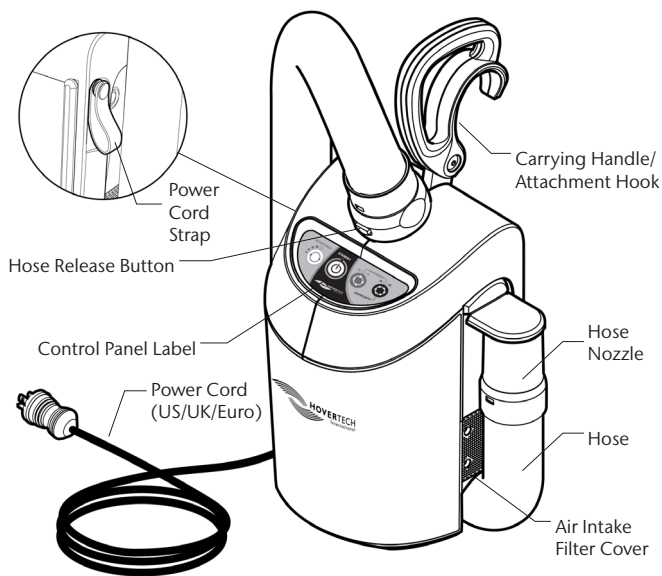

WARNING: The HT-Air is not compatible with DC power supplies.
The HT-Air is not for use with the Hoverjack Battery Cart.

HT-Air® Keypad Functions

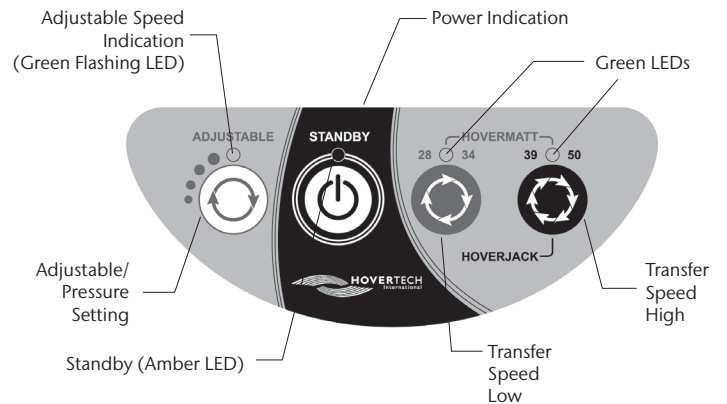

The ADJUSTABLE keypad function has four different settings. Each press of the button increases the air pressure and rate of inflation. The Green Flashing LED will indicate the inflation speed by the number of flashes (i.e. two flashes equals the second inflation speed).



STANDBY: Used to stop inflation/air flow (Amber LED indicates STANDBY mode).

Instructions for Use

To Laterally Turn a Patient:

1. Ensure the red deflation valves are capped tightly before beginning inflation.
2. Place the Q2Roller on the bed mattress and center using the Center Line Indicator.
3. Tuck the securing apron of the Q2Roller under the bed mattress to keep it in place during use.
4. Center an absorbent pad on top of the Q2Roller (if the patient is already on the mattress, place the device and pad under the patient using a log-rolling technique).
5. With the patient centered on the device, ensure all bed/stretchers rails are up, or follow your facility's protocol.
6. Locate the inflation valve on the OPPOSITE side of the turn direction, and place the air supply hose over the valve. Press the ADJUSTABLE button on the HT-Air keypad to initiate air flow. If necessary, press this button up to 3 more times to increase air flow and rate of inflation.
7. When the turn is complete, remove the hose from the valve and press the STANDBY button on the air supply keypad to stop air flow. Chamber may not need to be fully inflated to achieve desired position.
8. To quickly deflate the Q2Roller, remove the red cap from the deflation valve. To slowly deflate or adjust the chamber, depress the center of the black inflation valve.

To Use as a Post-Prone Positioning Device:

1. Follow steps 1-4 above.
2. Prepare the patient for proning.
3. Prone patient onto the Q2Roller and ensure patient is centered.
4. Inflate and adjust chambers to fully offload pressure and place patient in proper prone position per your facility's protocol.



Do not use Q2Roller to prone. Use other measures to prone patient - For Post-Prone Positioning only.

To Sit a Patient Up on the Edge of the Bed:

Use a minimum of 2 caregivers. Use clinical judgement to determine if the patient is capable and ready to sit on the edge of the bed.

1. Follow steps 1-4 to laterally turning a patient.
2. With 1 caregiver on each side of the patient, fully inflate the chamber on the opposite side of where patient will ultimately be seated to completely turn the patient. Bend the patient's knees for comfort.
3. Raise the head of the bed as high as it will go then lower the bed so that the patient's feet will touch the floor when seated at the edge of the bed. Both caregivers can now attend to the patient on the same side as the fully inflated chamber is supporting the patient's back.
4. One caregiver will assist with sliding the patient's legs to the side of the bed while the other is fully inflating the deflated chamber to ensure the patient does not slide off the edge of the bed.
5. Once the chamber is fully inflated and the patient is sitting upright, stabilize the patient and then open the deflate valve. As the chamber deflates, ensure the patient is still sitting upright. The leg caregiver will continue to guide the legs to the edge of the bed as the chamber is fully deflated. Stabilize the patient.

Product Specifications/Required Accessories

Q2ROLLER® LATERAL TURNING DEVICE

Material:	Polyurethane Laminated Non-Woven, Nylon Non-Woven
Construction:	RF Welded
Width:	Chambers: 46.5" (118 cm) (combined chamber width) Apron: 72" (183cm)
Length:	Chambers: 44.5" (113 cm) Apron: 44" (112 cm)

HoverCover Half Size

Material:	Non-Woven, Breathable PE
Construction:	Thermal Bonding
Width:	40" (101 cm)
Length:	40" (101 cm)

Model #: HTR-CHX-CASE (50 ea.)

Model #: HTR-200



LIMIT 600 LBS/ 272 KG



Single-Patient Use



Latex Free

REQUIRED ACCESSORY:

Model #: HTAIR1200 (North American Version) – 120V~, 60 Hz, 10A

Model #: HTAIR2300 (European Version) – 230V~, 50 Hz, 6A

Model #: HTAIR1000 (Japanese Version) – 100V~, 50/60 Hz, 12.5A

Model #: HTAIR2356 (Korean Version) – 230V~, 50/60 Hz, 6A

Cleaning & Maintenance

CLEANING

The Q2Roller is a single-patient use device. If the device becomes soiled, completely wipe it down using a germicidal cleaner (phenolic disinfectant, quaternary solution or other intermediate level disinfectant according to facility procedure) to disinfect. For best results, follow the cleaner manufacturer's recommended dwell time and instructions for use. Apply germicidal cleaner directly into hard-to-reach areas. Let air dry. Do not launder.

INFECTION CONTROL

The HoverCover Half Size is recommended to cover the Q2Roller to keep the device from getting soiled (also available for separate purchase). Other materials, such as pads or linens, may also be placed on top of the Q2Roller to keep it clean.

If the Q2Roller is used on an isolation patient, the hospital should employ the same protocols/procedures it utilizes for the bed mattress and/or for the linen in that patient room.

AIR SUPPLY CLEANING AND MAINTENANCE

See air supply manual for reference.

NOTE: CHECK YOUR LOCAL/STATE/FEDERAL/INTERNATIONAL GUIDELINES BEFORE DISPOSAL.

PREVENTIVE MAINTENANCE

Prior to use, a visual inspection should be performed on the Q2Roller to ensure that there is no visible damage that would render the Q2Roller unusable.

The Q2Roller should be periodically inspected to ensure the following:

- All inflation valves are self-sealing with no evidence of leakage
- There are no punctures or tears in the Q2Roller

If any damage is found that would cause the Q2Roller not to function as intended, the Q2Roller should be discarded.

Returns and Repairs

All products being returned to HoverTech must have a Returned Goods Authorization (RGA) number issued by the company. Please call (800) 471-2776 and ask for a member of the RGA Team who will issue you an RGA number. Any product returned without an RGA number will cause a delay in the repair time.

Returned products should be sent to:

HoverTech
Attn: RGA # _____
4482 Innovation Way
Allentown, PA 18109

For European companies, send returned products to:



Attn: RGA # _____
Kista Science Tower
SE-164 51 Kista, Sweden

For product warranties, visit our website:
<https://hovermatt.com/standard-product-warranty/>



HoverTech

4482 Innovation Way
Allentown, PA 18109

www.HoverMatt.com
Info@HoverMatt.com

These products comply with the standards applicable for Class 1 products in the Medical Device Regulation (EU) 2017/745 on medical devices.



CEpartner4U, ESDOORNLAAN 13,
3951DB MAARN, THE NETHERLANDS.

www.cepartner4u.com



Etac Ltd.

Unit 60, Hartlebury Trading Estate,
Hartlebury, Kidderminster,
Worcestershire, DY10 4JB
+44 121 561 2222

www.etac.com/uk



TapMed Swiss AG

Gumprechtstrasse 33
CH-6376 Emmetten
CHRN-AR-20003070

www.tapmed-swiss.ch

In case an adverse event in relation to the device, incidents should be reported to our authorized representative. Our authorized representative will forward information to the manufacturer.



4482 Innovation Way
Allentown, PA 18109

800.471.2776
Fax 610.694.9601

www.HoverMatt.com
Info@HoverMatt.com