




# POSITIONING WEDGES

When used together, the HoverTech Positioning Wedges and HoverMatt® Air Transfer System provide the optimal solution for pressure injury prevention and caregiver safety. With high density foam and a cover specially designed to stay in place under the patient, the positioning wedges effectively offload the sacrum, while the HoverMatt can be used to safely position the patient, eliminating skin shear and reducing the force required to transfer and position patients.

# POSITIONING WEDGES

HoverTech Positioning Wedges come in Single-Patient Use and Reusable styles, which both offer:

- High density foam to support the patient, while minimizing additional pressure from compression over time
- Anti-microbial, hypo-allergenic and latex-free foam
- 30 degree positioning complies with NPIAP Guidelines for the prevention and treatment of pressure injuries
- No-slip cover to secure and stabilize the device under the patient

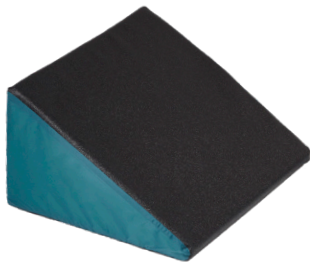

## SINGLE-PATIENT USE

- Polyurethane and nylon cover simplifies cleaning



## REUSABLE

- Dartex® cover with heat-sealed seams for added infection control
- Contour cut foam permits air circulation between the wedge and the patient
- Bariatric version comes with convex surfaces on the patient and caregiver sides for more support, and an extra-firm core topped with extra-plush memory foam for enhanced comfort and pressure redistribution
- 3 year warranty for standard version, 1 year for bariatric

Additional design features of Reusable Wedges include:



*Dartex® cover with heat-sealed seams*




*Waterfall Flap covers the top half of zipper enclosure*



*Raised border behind the bottom zipper prevents fluid infiltration*



*Contour cut foam for standard wedges*

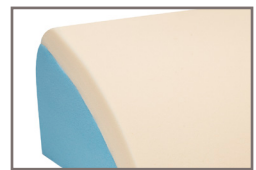

*Extra-firm core topped with extra-plush memory foam for bariatric wedge*

| STYLE                | MODEL NUMBER | ANGLE | SIZE              | QUANTITY |
|----------------------|--------------|-------|-------------------|----------|
| Single-Patient Use   | FPW-100      | 30°   | 15"L x 11"W x 7"H | Box - 2  |
| Reusable             | FPW-R-15S    | 30°   | 15"L x 11"W x 7"H | Each - 1 |
| Reusable             | FPW-R-20S    | 30°   | 20"L x 11"W x 7"H | Each - 1 |
| Reusable (Bariatric) | FPW-RB-26S   | 30°   | 26"L x 12"W x 8"H | Each - 1 |

4482 Innovation Way  
Allentown, PA 18109

800.471.2776

Fax 610.694.9601

[www.HoverMatt.com](http://www.HoverMatt.com)

[Info@HoverMatt.com](mailto:Info@HoverMatt.com)

BRO\_FPW\_090622

